



Sharpsmart Features



WIDE OPENING

One-hand, one-step deposition of sharps. The optimal access geometry of the lid opening, designed with a wide aperture for large and awkward sharps, minimises disposal-related injuries and restricts hand access to the contained waste.




NO LEAKS

A fixed liquid seal around the rim of the container lid eliminates leakage of bodily fluids during transport. Container leakage is a common hazard with sharps containers that are not equipped with a leakproof seal. Neoprene gasket in lid eliminates opportunity for blood or body fluid exposure



CANNOT BE REOPENED

Side locks are tamper proof when engaged. This safety feature prevents unauthorised access to the container contents. Temporary closure allows safe and easy transport within healthcare facility

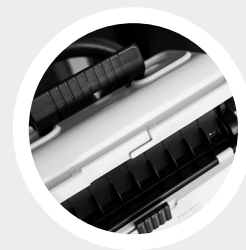

MAXIMUM STRENGTH HARDENED PLASTIC

The Medical grade hardened plastic casing of the Sharpsmart container is completely impenetrable by the contained sharps. The Sharpsmart container has been test-proven to an amazing 500 cycles proving its durability, sterility and re-usability.



ELIMINATES OVERFILLING

Once the container is full, the in-built safety tray is activated into an upright closed position to restrict access and overfilling. Passive overflow protection prevents staff from overfilling container and eliminates risk of overfilling injuries

Sharpsmart Sharps Containers

REGULAR



S14

14 Litre Regular Sharpsmart. Ideal for Patient Room use and lighter sharps volume areas

Fill Capacity: 6.5 Litres

H 285mm x W 398mm x D 210mm



S22

22 Litre Regular Sharpsmart. Best suited to moderate volume areas

Fill Capacity: 14.5 Litres

H 446mm x W 398mm x D 210mm



S32

32 Litre Regular Sharpsmart. Developed for high volume areas and for use with Cartsmart 1 & floor stand accessories

Fill Capacity: 23.5 Litres

H 618mm x W 398mm x D 210mm

ACCESS PLUS



S14 AP

14 Litre Access Plus Sharpsmart. Designed for secure areas where unrestricted deposit of lighter sharps volume is required

Fill Capacity: 6.5 Litres

H 285mm x W 398mm x D 210mm



S32 AP

32 Litre Access Plus Sharpsmart. Designed for secure areas where unrestricted deposit of long or bulky sharps is required

Fill Capacity: 32 Litres

H 618mm x W 398mm x D 210mm



S64 AP

64 Litre Access Plus Sharpsmart. Designed for secure areas generating large volumes of long or bulky sharps

Fill Capacity: 47 Litres

H 726mm x W 398mm x D 354mm